

OLIGONUCLEOTIDE FOR GENOTYPING *MYCOPLASMA* AND ITS
RELATED STRAINS, MICROARRAY COMPRISING THE
OLIGONUCLEOTIDE, AND METHOD FOR DETECTING STRAINS
USING THE MICROARRAY

5

Technical Field

The present invention relates to a method for detecting
Mycoplasma and its related strains which are a source of contamination of
cell lines and biological products and human pathogens. More particularly,
10 the present invention relates to genus-specific and species-specific
oligonucleotides for genotyping *Mycoplasma*, *Acholeplasm* and
Ureaplasma strains, a microarray comprising the oligonucleotides, and a
method for detecting strains using the microarray

15 **Background Art**

Mycoplasma is a prokaryote pertaining to Mollicute family without
cell wall, which was known as a hospital acquired pathogen causing
pneumonia via infection of genital and respiratory organs of human as
well as livestock such as pig and cow. Recently, *Mycoplasma* is more
20 seriously understood as a major contaminant of cell culture and cell line

Especially, as the development and production of biological
products for protecting and treating human diseases increases, the
contamination of various pathogens provided by microorganism or
clinical sample in the process of production became a serious problem.
25 Examples of the biological products are an oncolytic virus, vaccine, a
gene therapy vector and a recombinant protein. They have been found
to be contaminated by bacteria, fungus, virus, *Mycoplasma* and its
related strains (Doblhoff-Dier et al., 2001). The reason of the
contamination is an organism contaminated in media components or
30 experimental instruments and cross-contamination of microorganism and
virus in air (Jung et al., 2003). Also, the contamination can be occurred

by a cross-contamination of already-infected WCB (Working Cell Bank) which is used for mass production of biological products (Wisher et al., 2002).

It is reported that, among these contamination sources, about 15~35% of cell culture or cell line is infected by *Mycoplasma* and its related strains (Hopert et al., 1993). This also makes experimental results incredible because it can change characteristics of cells such as abnormal synthesis of DNA, RNA and protein by binding to host cell wall (Kong et al., 2001). As gene therapy and cell therapy are getting into the spotlight recently, an assay for infection of stem cell and cord blood by *Mycoplasma* and its related strain became more important. Therefore, for the credible and reproducible experimental results and the quality control of commercialized biological products, it is essential to detect an infection with *Mytco*plasma and its related strains.

Under this situation, Europe community make it a rule that, for credibility of safety and quality of food and drug, GMP (Good Manufacturing Practice) and QC (Quality Control) should be submitted and cell banks such as MCB (Master Cell Bank) and WCB should be subjected to an assay for detection of virus, fungus and bacteria such as *Mycoplasma* (Doblhoff-Dier et al., 2001).

About 100 kinds of bacteria pertaining to Mollicute family without cell wall have been found so far, including *Acholeplasma*, *Enteroplasma*, *Mesoplasma*, *Mycoplasma*, *Ureaplasma* and *Spiroplasma*. Among them, about 20 kinds of *Mycoplasma*, *Acholeplasma* and *Ureaplasma* are major contamination source of cell culture. These are referred to as "*Mycoplasma* and its related strains" in this specification. About 95% of the contaminants are covered by *M. arginini*, *M. fermentans*, *M. orale*, *M. hyorhina*, *M. hominis*, *M. salivarium*, *M. pirum*, *A. laidlawii* (Dorigo-zetsma et al., 1997). However, *Mycoplasma* is difficult to be cultured in extracellular media and turbidity is rare in the culture. Therefore, there has been a need to the rapid and accurate genotypic detection method

which can trace a contamination source of *Mycoplasma* and its related strains.

Conventional *Mycoplasma* detection methods are the culturing method, the DNA fluorochole stain method, the immunofluorescence method, and the polymerase chain reaction (PCR) method (Dorigo-
5 zetsma et al., 1997). However, the culturing method has a drawback that extracellular culturing is difficult, preparing its media is complex by adding supplements such as serum and culturing time is too long, about 4 days ~ 3 weeks according to the kinds of strains (Jensen et al., 2003).

10 The DNA fluorochole stain method such as Hoechst 33258 stain has a drawback that culturing condition is too difficult to match and subjective inspectors can make a misjudgment (Chen et al., 1997). The immunofluorescence method such as ELISA has a drawback that bacteria having similar antigen with *Mycoplasma* such as *Streptococcus*
15 *milleri* group and *Staphylococcus aureus* may raise a false positive signal due to of low specificity (Hopert et al., 1993). The PCR method makes use of 16S/23S intergenic spacer region (ITS) and a gene coding 169 kDa of P1 cyadhesion proteine which represent variety of *Mycoplasma* (Uphoff et al., 2002). The P1 gene, a surface antigen
20 gene, has several subtypes representing diversity and has been used as a target gene for serological detection using immune reaction and genotypic detection using restriction fragment length polymorphism (RFLP) to identify *Mycoplasma* (Campo et al., 1998). However, most of conventional PCR methods use a primer designed based on 16S rRNA
25 which is a common sequence of prokaryotes, and second PCR or nested PCR having high sensitivity can make a cross-contamination of *Mycoplasma* dispersed in air and an amplification of a bacteria similar with *Mycoplasma* in classification (Uphoff et al., 2002).

To overcome the above limitations of the conventional detection
30 methods, a genotypic detection method using probes have been developed recently, which make it possible to analyze many kinds of

genes in a short time using DNA hybridization principle based on gene sequencing and detect specifically a single base change using a proper hybridization condition between specific probe and target DNA.

5 The present inventors developed ITS-derived oligonucleotides capable of detecting *Mycoplasma* and its related strains, which are important in genotypic detection, and a microarray comprising the oligonucleotides as a probe for detecting *Mycoplasma* and its related strains.

10 **Disclosure of the Invention**

It is a first object of the present invention to provide oligonucleotides for detecting *Mycoplasma* and its related strains designed based on their ITS base sequences.

15 It is another object of the present invention to provide novel ITS sequences of *Mycoplasma bovis*, *Mycoplasma cloacale*, *Mycoplasma falconis*, *Mycoplasma faucium*, *Mycoplasma spermatophilum* and *Mycoplasma synoviae*, which is useful for detecting *Mycoplasma* and its related strains.

20 It is another object of the present invention to provide a microarray comprising genus-specific and species-specific oligonucleotides for detecting *Mycoplasma* and its related strains as probes.

It is another object of the present invention to provide a method for detecting *Mycoplasma* and its related strains using the microarray.

25 It is another object of the present invention to provide a kit for diagnosing *Mycoplasma* and its related species infection individually or simultaneously, comprising genus-specific and species-specific oligonucleotides for genotyping *Acholeplasma* and related strains.

30 According to an aspect of the present invention, there is provided a purified ITS (internal transcribed spacer) target DNA for genotyping

Mycoplasma strains, comprising any one sequence selected from SEQ ID Nos. 1 to 6.

SIQ ID Nos. 1 to 6 are base sequences of ITS (internal transcribed spacer) of *Mycoplasma bovis*, *Mycoplasma cloacale*,
5 *Mycoplasma falconis*, *Mycoplasma faucium*, *Mycoplasma spermatophilum* and *Mycoplasma synoviae*, which was newly obtained by base sequencing analysis.

The ITS target DNA of the present invention can be used indirectly for designing probes or primers used for genotyping
10 *Mycoplasma* strains or directly for genotyping *Mycoplasma* strains via PCR amplification.

According to another aspect of the present invention, there is provided an oligonucleotide for genus-specific genotyping of *Mycoplasma* and *Ureaplasma* strains, comprising any one sequence
15 selected from SEQ ID Nos. 7 to 21 or its complementary sequence.

According to another aspect of the present invention, there is provided an oligonucleotide for genus-specific genotyping of *Acholeplasma* strains, comprising any one sequence selected from SEQ ID Nos. 22 to 27 or its complementary sequence.

20 According to another aspect of the present invention, there is provided an oligonucleotide for species-specific genotyping of *Mycoplasma* and *Ureaplasma* strains, comprising any one sequence selected from SEQ ID Nos. 28 to 127 or its complementary sequence.

According to another aspect of the present invention, there is provided an oligonucleotide for species-specific genotyping of *Acholeplasma* strains, comprising any one sequence selected from SEQ
25 ID Nos. 128 to 133 or its complementary sequence.

The oligonucleotides according to the present invention are designed based on multiple sequence alignment of ITS (internal
30 transcribed spacer) sequences, which are present between 16S rRNA and 23S rRNA of *Mycoplasma* and its related species. The

oligonucleotides can be used as primers for PCR amplification in order to genotype *Mycoplasma* and its related species or as probes for hybridization reaction in order to genotype *Mycoplasma* and its related species.

5 According to another aspect of the present invention, there is provided a microarray comprising more than one oligonucleotides selected from genus-specific and species-specific oligonucleotides for genotyping *Mycoplasma*, *Acholeplasma* and *Ureaplasma* strains according to any one from claims 2 to 5 as probes attached on a support.

10 In the microarray according to the present invention, the probes may be any materials having base sequence, preferably any one selected from a group consisting of DNA (Deoxyribose Nucleic acid), RNA (Ribose Nucleic Acid), and nucleic acid analogues such as PNA (Peptide Nucleic Acid), LNA (Locked Nucleic Acid) and HNA (Hexitol
15 Nucleic Acid).

 In the microarray according to the present invention, the support may be any materials to which the probes can be attached, preferably any one selected from a group consisting of slide glass, plastic, membrane, semiconductive chip, silicon and gel. The microarray
20 according to the present invention can be manufactured using conventional method such as pin microarray, ink jet, photolithography or electric array method.

 The microarray according to the present invention can be used for simultaneously genotyping various *Mycoplasma* and its related
25 species which are known as a major contaminant of biological drug and cell line as well as a human pathogen from one sample, as the microarray comprises genus-specific and species-specific oligonucleotides for genotyping *Mycoplasma* and its related species as a set attached a support

30 According to another aspect of the present invention, there is provided a method for detecting *Mycoplasma*, *Acholeplasma* and

Ureaplasma strains, comprising the following steps:

- a) extracting nucleic acids from a sample;
- b) amplifying target DNA among the extracted nucleic acids;
- c) hybridizing the amplified target DNA with probes of the
- 5 microarray according to the above present invention; and
- d) detecting signals generated from the hybridization reaction.

In the detection method according to the present invention, the sample may be biological drug, cell line, or human tissues or serum. The purifying step can be performed using conventional DNA or RNA
10 purification method or kit. The signal detecting step can be performed using a conventional fluorescence scanner after binding conventional fluorescent dyes such as Cy5 or Cy3.

According to another aspect of the present invention, there is provided a kit for diagnosing *Mycoplasma* and its related species
15 infection, comprising more than one oligonucleotide selected from genus-specific and species-specific oligonucleotides for genotyping *Acholeplasma*, *Mycoplasma* and *Ureaplasma* strains according to the above present invention.

In the kit according to the present invention, the oligonucleotides
20 are used as probes for hybridizing with target sample and may be contained in a proper vessel. The probes may be labeled with a radioactive or non-radioactive labeling agent, the latter comprises conventional biotin, Dig(digoxigenin), FRET(fluorescence resonance energy transfer) or fluorescent dye (Cy5 or Cy3). Further, the
25 oligonucleotides can be used as primers for PCR amplification. In this case, the kit may contain DNA polymerase, 4 dNTPs and PCR buffer for PCR reaction. In addition, the oligonucleotides can be attached to a microarray as probes. In this case, the kit may contain hybridization reaction buffer, PCR kit containing primers for amplifying a target gene,
30 washing solution for the unhybridized DNA, dyes, washing solution for unbound dyes and manual sheet for the microarray.

Hereafter, the present invention will be described in more detail.

The present invention provides a method for detecting or genotyping *Mycoplasma* and its related strains which is a major contamination source of cell lines and biological products and a human pathogen, comprising the following steps:

a) if necessary, extracting nucleic acids from a sample such as cell lines, biological products or human tissue or serum;

b) if necessary, amplifying target DNA of *Acholeplasma*, *Mycoplasma* and *Ureaplasma* strains among the extracted nucleic acids using more than one proper primers;

c) hybridizing the amplified target DNA with probes having a sense or antisense or complementary sequences of genus-specific and species-specific oligonucleotides of *Acholeplasma*, *Mycoplasma* and *Ureaplasma* strains disclosed in Tables 2 and 3; and

d) detecting signals generated from the hybridization reaction.

From the detected signals in the step d), the existence of *Mycoplasma* and its related strains in the sample can be predicted.

The present inventors carried out a sequence analysis of ITS regions of many *Acholeplasma*, *Mycoplasma* and *Ureaplasma* strains to obtain genus-specific and species-specific oligonucleotides for detecting *Mycoplasma* and its related stains which can be a basis of developing a specific and sensitive hybridization assay. Also, the present inventors newly analyzed ITS sequences of newly found 6 *Mycoplasma* strains, which makes it possible to design probes capable of detecting more various *Mycoplasma* and its related strains.

Table 1 discloses ITS sequences of newly analyzed 6 strains among target sequences for detecting *Mycoplasma* strains, which correspond to SEQ ID Nos. 1 to 6. In the present invention, the probes for detecting *Mycoplasma* strains were designed based on the multiple alignment of ITS sequences of *Mycoplasma*.

FIGS. 1 and 2 show multiple sequence alignments of ITS regions

of *Mycoplasma*, *Acholeplasma* and *Ureaplasma* for selecting genus-specific and species-specific probes of *Mycoplasma* and its related strains. Genus-specific oligonucleotides of *Mycoplasma* and *Ureaplasma* were designed from conservative sequence region indicated by a box in FIGS. 1a to 1f. Species-specific oligonucleotides of *Mycoplasma* and *Ureaplasma* were designed from polymorphic sequence region outside the box in FIGS. 1a to 1f. Genus-specific oligonucleotides of *Acholeplasma* were designed from conservative sequence region indicated by a box in FIGS. 2a to 2c. Species-specific oligonucleotides of *Acholeplasma* were designed from polymorphic sequence region outside the box in FIGS. 2a to 2c.

In step b) of the present invention, the target DNA of *Acholeplasma*, *Mycoplasma* and *Ureaplasma* strains were amplified using more than one pair of proper primers. FIG. 3 shows PCR amplification of ITS target sequences of *Mycoplasma* and its related strains using a primer pair, MP16SF-2 and MP23SR-2. In FIG. 3, 1 is a PCR product of *M. arginini*, 2 is a PCR product of *M. arthritidis*, 3 is a PCR product of *M. fermentans*, 4 is a PCR product of *M. hominis*, 5 is a PCR product of *M. hyorhinitis*, 6 is a PCR product of *M. neurolyticum*, 7 is a PCR product of *M. opalescens*, 8 is a PCR product of *M. orale*, 9 is a PCR product of *M. pirum*, 10 is a PCR product of *M. penetrans*, 11 is a PCR product of *M. pulmonis*, 12 is a PCR product of *M. salivarium*, 13 is a PCR product of *M. cloacale*, 14 is a PCR product of *M. falconis*, 15 is a PCR product of *M. faucium*, 16 is a PCR product of *M. hyosynoviae*, 17 is a PCR product of *M. muris*, 18 is a PCR product of *M. primate*, 19 is a PCR product of *M. spermatophilum*, 20 is a PCR product of *M. synoviae*, 21 is a PCR product of *M. pneumoniae*, 22 is a PCR product of *M. genitalium*, 23 is a PCR product of *M. bovis*, 24 is a PCR product of *U. urealyticum*, 25 is a PCR product of *A. laidlawii*.

In step c) of the present invention, the amplified target DNA were hybridized with probes for detecting *Mycoplasma* and its related strains.

Preferably, the probes may be a combination of more than one probes capable of simultaneously detecting many *Mycoplasma* and its related strains from a single sample. Practically, the probes are optimized to simultaneously hybridize with multiple target DNAs of *Mycoplasma* and its related strains under the same hybridization and washing conditions.

The present invention provides a microarray comprising a set of probes for detecting *Mycoplasma* and its related strains, which can simultaneously detect many *Mycoplasma* and its related strains from a single sample with a single experiment.

In the present invention, the term 'probe' means a single-stranded oligonucleotide having a sequence complementary to target DNA of *Mycoplasma*, *Acholeplasma* and *Ureaplasma*. The probe may have a sense, antisense or complementary sequence of SEQ ID Nos. disclosed in this specification as long as it can hybridize with one of double strands of target DNA. The oligonucleotide may be ribonucleotide (RNA), deoxynucleotide (DNA), peptide nucleic acid (PNA) or locked nucleic acid (LNA), and contain modified nucleotides such as Inosine only if it does not change their hybridization characteristics. Preferably, the genus-specific oligonucleotides for detecting *Mycoplasma*, *Acholeplasma* and *Ureaplasma* may have a base sequence of SEQ ID Nos. 7 to 27. Preferably, the species-specific oligonucleotides for detecting *Mycoplasma*, *Acholeplasma* and *Ureaplasma* may have a base sequence of SEQ ID Nos. 28 to 133.

FIG. 4 shows a microarray comprising probes for detecting genotypes of *Mycoplasma* and its related strains as a set on a support. In FIG. 4, each species name and SEQ IN Nos. are described which correspond to individual probes. The terms 'MP-C' and 'AP-C' mean *Mycoplasma* and *Ureaplasma* genus and *Acholeplasma* genus. FIG. 4 is no more than an example of probe compartment of the present invention, so compartment and layout of each probe can be varied.

In the present invention, newly analyzed ITS sequences of 6

Mycoplasma strains as a target DNA for detecting *Mycoplasma* and its related strains are as shown in Table 1. The genus-specific oligonucleotides for detecting *Mycoplasma*, *Acholeplasma* and *Ureaplasma* used in the present invention are as shown in Table 2. The species-specific oligonucleotides for detecting *Mycoplasma*, *Acholeplasma* and *Ureaplasma* used in the present invention are as shown in Table 3.

【Table 1】

Species	Sequence (5' → 3')	SEQ ID NO.
<i>M.bovis</i>	TTC TACGGAGTACACTTGTCTTTTATCACTATAAAAAAAGAC TTATAACCAAAAT TAC TAGACC TATATTTATTTATAAACGTCATGGCTTTTATTAATAGG TCAAAAGCTA TATATCTAGTTTGTAGAGAACATTC TCATATGTTC TTGTAAAAC TGAA TAG TAAA ATA TTTTTCGATATTTACAACGACATCAAAAA TCAAA TTAATGG TTAATTTG TTTTG ATTCA TCG AG TAAGTCATATTTAATA TGTTCATTGAAA TGTCTTAAAA TACACATC TAAAC TAACAAC AATAGGAAAA TACTAC TTTTAAA TAAGGAAGAG TTTTGTGTGG ATGC	1
<i>M.cloacae</i>	CTTC TACGGAGTACAATTC TCAC TG TTATGGAA TTAATTTGTATCCAGTTTGTAGA GAAC TTTCTCTCAATTTTGTTC TTGTAAAAAC TGAA TATAGACATTG AAA TCAATAAA TTAATATTTCAAA TGT TTAGATCAACCTATAGAAATATCAAGACATATACAAAAATA GGTCATAC TTATATTTATAAATAC T	2
<i>M.falconis</i>	CTTTC TACGGAGTACAAC TTC TG TTATGGAA TAA TATTGTATCCAG TTTGTAGAGT ACTAATCTCTCTTTTGTTC TTGTAAAAAC TGAA TATCGACATTGAAAAATTA TTAAT AATATTTC AAAG TTTAGATCAACCTATAGAA TACAAAAATA TAGACAACAA TAGGT CATACAACAAACA TAACAAAACAAC T	3
<i>M.faucium</i>	GAA TGG TGGCTTCGAGACTAAAAG TTA TG GAAAAACA TCG TATCCAG TTTGTAGA GAAC TAAAC TTC TC TTTTGTTC TTGTAAAAAC TGAA TATAGACATTGAAAAATTA AAA TTAATA TTTCAAAG TTTAGATCAACCTATAG AATACAAAA TCAATACAATAGG TCAATAC TATACAATTGCA TAACAAAAA TACTATTAACAAGATAAGAG TTTTGT GTGGATGCAATTGTAT	4
<i>M.sppermatophilum</i>	GTGGGGATGGATCACC TCC TTTCTACG GAGTACAACATACATTCAAA TTTGACT GAA TG TTA TTAACC TTA TTTTTCAC TAGGCC TTTTAA TATATTTGT TATG TGACT TTTATGGCTTAAAAGTCTTATATCTAG TTTTGTAGAGGACATCC TCCTAA TTGTTC TTGAAAAC TGAATAG TAAATTTTGTATATTTACAACGACATC TAAATAA TTGAA TT AAGTCAATTTGT TTAGATTTCA TCGAGATAG TCA TTTTAAAAAATGA TTCA TTGAA ATGTCTTAAAATACACA TCAAAACAAACAATC TATACAATAGGAA TTTATATACT	5
<i>M.synoviae</i>	TCC TTACG GAGTACATTAATTTTACAAAAGGCATTTTATTAAC TG AAAGC TTTTAG AGAAAA TTCTAAAAGCGG TTGTGTATCGC TTTTGTGCTTGGGC TATTGTATTTA GTTTGTAGAGAACAAACC TC TC TTA AAAA TTGTTC TTGTAAAAC TAAA TAG TAA TAAA GATATTACAACGACATCAAAAA TATAA TTAATTAAGGTAA TTTGT TTTTGA TACCG AG TTTAAA TTA TGT AATAA TAA TTTATTA AAA TGTCTTGTAA TACATCA TAACAAATA TAACAA TAGGACATATGTATAC TAAC TTTTAAAAAGT	6

【Table 2】

Genus	Probe	Sequence	SEQ ID NO.
<i>Mycoplasma</i>	MP-CP1	TTCTTTGAAAAC TGA	7
	MP-CP2	FWTC TTTVAAAAC TTTATWN	8
	MP-CA1 MP-CA2	MWTYG TTTCCAG TTTGAGAG TTTAGATCAACCTA TAGAATA	9
			10
	MP-CB1 MP-CB2 MP-CB3 MP-CB4	RTATYTAG TTTGAGAGRRCA WWIRATTYATTAAATGCTT GGKYAATTGTTTGGAT RATATTTACAMCGMCAYC	11
			12
			13
			14
	MP-CC1 MP-CC2	CCTCCTTCTATCGGAGTAMA CGGATTCTATTTAGTTTGAG	15
			16
	MP-CD1 MP-CD2 MP-CD3	TAAAA TAGATACCTTAAKATA GTATYYAGTTTGGAAAG CTTGCCAAWTAGWTWT	17
			18
			19
	MP-CE1 MP-CE2	AWACHACAACTTTCTAGTTC AATAAGTTACTAAGGGCTTAT	20
			21
<i>Acholeplasma</i>	AP-CP1	TCA TCATATTCAGTTTIG	22
	AP-CA1	GGGCC TTTAGCTCAGYTGGTT	23
	AP-CA2	AGAGCFCWCGCYTGA TAAGCG	24
	AP-CA3	WGRGGTCGATGGTTCHAGTCC	25
	AP-CB1	TCATCATATTCAGTTTIGARRH	26
	AP-CB2	AGTC TTTGAAAAGTAGATAAA	27

【Table 3】

Species	Probe	Sequence	SEQ ID NO.
<i>M. arginini</i>	MP-arg1	AGATTATATCATACAATAGA	28
	MP-arg2	GAGTACATAAATGTTATGGAA	29
<i>M. arthritidis-faucium</i>	MP-arf1	TGAAGCCCGATGGTGCTTCG	30
	MP-arf2	TGAGAGAAC TAAACTTC TCTC	31
	MP-arf3	GAATACAAAA TCAATACAATA	32
<i>M. fermentans</i>	MP-fer1	ATGTACTATTAAC TTATTCAC	33
	MP-fer2	TACAAAAGAGTAC TTTTAAAA	34
	MP-fer3	TTTTATGGGTC TAAAGCTTT	35
	MP-fer4	GAACAATATTTTTTCTC TCA	36
	MP-fer5	ATAACAACTATAACAATAGG	37
<i>M. hominis</i>	MP-hom1	ATTATCTCTCGGTCTTT	38
	MP-hom2	ATATTTATATTTTATAAGACA	39
	MP-hom3	ATTGATATATTAATTAAATATT	40
<i>M. hyorthinis</i>	MP-hyo1	GAATAGCAAATAACAATATGATT	41
	MP-hyo2	CGGAGTACATTAGTCTTAATT	42
	MP-hyo3	TTACATAATCGATTCTGTCT	43
	MP-hyo4	AGCTTTAAGTCTCAATTATA	44
	MP-hyo5	TTCAATTTATTTATTTCAACG	45
	MP-hyo6	AACGATCTTTTTTATAAACC GA	46
	MP-hyo7	TTAAATTTCTAAATAGATTA	47
	MP-hyo8	AGATATTTATCTTTAGCAATA	48
<i>M. neurolyticum</i>	MP-neu1	GGTTATTATGGGCTTGCTA	49
	MP-neu2	GGTTATTTAAAAATCC TTTTA	50
	MP-neu3	TAATTTTTCTTTCTAAATIAA	51
<i>M. opalescens</i>	MP-opa1	CATCATAATGTAACCAATAC	52
	MP-opa2	ACAAAAATCATTA TTTTAAAT	53
	MP-opa3	TTTAATGATTATTAACCTTTT	54
	MP-opa4	TTATGTGCTTTG TTTTATGG	55
	MP-opa5	TATGGTC TACAAAGCTTATAT	56
	MP-opa6	GATAAAAAACAATCATAAATT	57
<i>M. orale</i>	MP-ora1	CATAAATAGTTAATGGCTCA	58
	MP-ora2	ATAGAGACAAATACAAAAACA	59
	MP-ora3	GGTCAAAAAATACTTATACGTA	60
<i>M. pirum</i>	MP-pir1	TAGTCTTTGGTGGAATAACA	61
	MP-pir2	CTTTATACACCTTATTACAAT	62
	MP-pir3	TAAAAATCCAATTTAAATGTTA	63
	MP-pir4	GCAAATTTGATGTCACATTT	64
	MP-pir5	AATTAATCTCTCCTATTACTT	65

	MP-pir6	TTAAAGTAGTAGAGATGGTTC	66
	MP-pir7	CAAATATCAAAATGCTAATGGA	67
	MP-pir8	ATGC TAA TGG ATATCAAAAAA	68
<i>M. penetrans</i>	MP-pen1	AAGAGTAAGTTC TAGGTCG	69
	MP-pen2	CATTAAAGCTAAGTAACAAAT	70
	MP-pen3	TCC TAAAC TGAAATTTATCT	71
	MP-pen4	TTATATAAGAGTAAGTTC TAG	72
	MP-pen5	ATTTTTC TC TCAAGATAGTTC	73
	MP-pen6	TC TAATCATACTTG TTATTTT	74
<i>M. pulmonis</i>	MP-pul1	AATTTTGTATCCGAGTCA TT	75
	MP-pul2	CATTTTCTATCAATAGTTAT	76
	MP-pul3	TATG TGTATC TTGCCAATTAG	77
	MP-pul4	TTCTATCTTTCAAACAAATA	78
	MP-pul5	TATAAATTAAATATGATAACG T	79
	MP-pul6	TCATCAAAATG TAAAAATTTT	80
	MP-pul7	AAAAA TAAAA TAGATACCTTA	81
	MP-pul8	AAATAAA TTCAACAATAGGA	82
<i>M. salivarium</i>	MP-sal1	TAATGGATTAAATTTTCGTG	83
	MP-sal2	TATCAAAATCAATATAATATTT	84
<i>M. cloacae</i>	MP-clo1	AGTACAAATCTCAC TGTATG	85
	MP-clo2	TAGAATATTCAGACATATAC	86
<i>M. falconis</i>	MP-fal1	GAGTACAAC TTCTGTTATG	87
	MP-fal2	AGAA TACAAAAATATAGACAA	88
	MP-fal3	ATTGAAAAATTAATTAAT	89
<i>M. hyosynoviae</i>	MP-hyos1	CTAGACTAAAGTTAA TGGTAC	90
	MP-hyos2	AA TTA TCAAATTAATTTCA	91
<i>M. muris</i>	MP-mur1	TATAGAAAACCCACATCA	92
	MP-mur2	TATAGAATA TTTTAAATTT	93
	MP-mur3	GATTA TTACACCA TATTAGAA	94
	MP-mur4	TCAATAAACCTAAATAAAAA	95
<i>M. primate</i>	MP-pri1	GTAGACATAACCCAGCTA	96
	MP-pri2	CAAACGTCTATCGCTTTT TAG	97
	MP-pri3	TCATGGGCTTTTAATAGGGTC	98
	MP-pri4	ACCCCAAC TCCCATCAAAAAAT	99
<i>M. spermatophilum</i>	MP-spe1	TTATCGAGATAGTCATTTTA	100
	MP-spe2	CAAACATACATTCAAATTTT	101
	MP-spe3	TTTGTACTGAA TGTATTAAC	102
	MP-spe4	TTTGTATG TGACTTTTATGG	103
	MP-spe5	AAAACAAACAATCTATACAAT	104

<i>M. synoviae</i>	MP-syn1	TTGGCTTGGGCTATTGTATT	105
	MP-syn2	GCGGTTGTGTATCGCTTTTT	106
	MP-syn3	ACCTCTCTAAAAATTGTTCTT	107
	MP-syn4	CCGAGTTTAAATTATTGAATA	108
	MP-syn5	CATCATACAACAATAACAATA	109
<i>M. pneumoniae</i>	MP-pne1	GTAATTAAACCCAAATCCC	110
	MP-pne2	ATCTTTAATAAGATAAATAC	111
	MP-pne3	CTAAACAAAACATCAAAATCC	112
	MP-pne4	AAAGAACATTTCGCTTCTTT	113
<i>M. genitalium</i>	MP-gen1	CACCCCTTAATTTTTTCGG	114
	MP-gen2	AATGGAGTTTATTATTTTATTAA	115
	MP-gen3	CCCAATCAATGTTTGGTCCTC	116
	MP-gen4	CAACTAACACACTTGGTCAGT	117
	MP-gen5	AGAATGTTTTTGAACAGTTC	118
	MP-gen6	TAGTCCAAAAATAAATACCA	119
<i>M. bovis</i>	MP-bov1	TATAACCAAAATTAAGACCTA	120
	MP-bov2	GTCATGGCTTTATTAAATAGG	121
<i>U. urealyticum</i>	UP-ure1	CATTAAAGTTGTCAGTGAA	122
	UP-ure2	TAATTTACGTACTAATAAGTG	123
	UP-ure3	TTTATTAAAAATCCATATGAAT	124
	UP-ure4	AAGCCACTTTTTTAAAAATTT	125
	UP-ure5	CCATAATAATTAAATTATTAT	126
	UP-ure6	ATTATCAACAAATCTTTCATA	127
<i>A. laidlawii</i>	AP-lai1	AACACTTAGCACAAAGATGAC	128
	AP-lai2	CTTCTAAGGAGAAAGGCATA	129
	AP-lai3	ATGACTACTAGTAAGTAGTAA	130
	AP-lai4	GTAGTAATAATTCCTAAATTT	131
	AP-lai5	TTAAAGTAATTTAAGTGTTTC	132
	AP-lai6	TAAATGATGCTGAAAAGAAA	133

* Mixed Base \Rightarrow Code Name

M : A + C, W : A + T, Y : C + T, R : A + G
K : G + T, V : G + A + C, N : A + G + C + T

5

Brief Description of the Drawings

FIGS. 1a to 1f show multiple sequence alignments of each ITS region of *Mycoplasma* and *Ureaplasma* for selecting genus-specific probes.

FIGS. 2a to 2c show multiple sequence alignments of each ITS region of *Acholeplasma* for selecting genus-specific probes.

FIG. 3 shows a result of PCR amplification using primer pairs which can amplify ITS target sequences of many *Mycoplasma* and its related strains

FIG. 4 shows a microarray comprising probes for detecting genotypes of *Mycoplasma* and its related strains as a set on a support.

FIGS. 5a to 5k show results of image analysis of specific hybridization reaction of each probes for detecting genotypes of *Mycoplasma* and its related strains and results of numerical analysis calculated from pixel intensity.

5

Best mode for carrying out the Invention

The present invention will be described in greater detail by means of following examples. The following examples are for illustrative purpose and are not intended to limit the scope of the invention.

10

Example 1: Incubation of *Mycoplasma* and its related strains and Isolation of Genomic DNA

Total 25 kinds of strains, including 1 kind of *Acholeplasma*, 23 kinds of *Mycoplasma*, and 1 kind of *Ureaplasma* were obtained from the American Type Culture Collection (ATCC). The strains were cultured in each culturing media under each culturing conditions according to manual provided by ATCC. From the cultured media, strain colonies were obtained with a white gold ear and input in 1.5ml tube, 100μl of InstaGene matrix (Bio-Rad, USA) was added thereto and suspended, and reaction was performed at 56°C for 30 minutes in constant temperature bath. And then, the reactant was shook for 10 seconds, heated at 100°C for 8 min, shook again for 10 sec, centrifuged at 12,000 rpm for 3 min, transferred to new tube, and freeze-stored at -20°C. The product was used as template DNA of PCR reaction.

25

The strains used were as followed:

Acholeplasma laidlawii (ATCC 25937)

Mycoplasma arginini (ATCC 23838)

Mycoplasma arthritidis (ATCC 19611)

Mycoplasma bovis (ATCC 27368)

30

Mycoplasma cloacale (ATCC 35276)

Mycoplasma falconis (ATCC 51372)
Mycoplasma faucium (ATCC 25293)
Mycoplasma fermentans (ATCC 19989)
Mycoplasma genitalium (ATCC 33530)
5 *Mycoplasma hominis* (ATCC 23114)
Mycoplasma hyorhinis (ATCC 17981)
Mycoplasma hyosynoviae (ATCC 25591)
Mycoplasma muris (ATCC 33757)
Mycoplasma neurolyticum (ATCC 19988)
10 *Mycoplasma opalescens* (ATCC 27921)
Mycoplasma orale (ATCC 23714)
Mycoplasma penetrans (ATCC 55252)
Mycoplasma pirum (ATCC 25960)
Mycoplasma pneumoniae (ATCC 15531)
15 *Mycoplasma primum* (ATCC 15497)
Mycoplasma pulmonis (ATCC 14267)
Mycoplasma salivarium (ATCC 23064)
Mycoplasma spermatophilum (ATCC 49695)
Mycoplasma synoviae (ATCC 25204)
20 *Ureaplasma urealyticum* (ATCC 27618)

Example 2: Preparation of probes for detection of *Mycoplasma* and its related strains

The probes used for detection of *Mycoplasma* and its related
25 strains were selected based on a result of multiple alignment of ITS
sequences of *Mycoplasma*, *Acholeplasma* and *Ureaplasma*. Among
Mycoplasma and its related species, 16S rRNA sequences has high
similarity of 74~97%, whereas ITS sequences has lower similarity of
25.4~78.8% except for between *M. salivarium* and *M. hyosynoviae*, and
30 *M. hominis* and *M. falconis*. In other words, ITS contains a region more
polymorphic than 16S rRNA which is useful for designing probes for

detection of *Mycoplasma* and its related strains. However, to complement specificity between *M.salivarium* and *M. hyosynoviae*, and *M. hominis* and *M. falconis* having a high ITS similarity, more restrictive and strict probes were designed.

5 In the present invention, the oligonucleotide probes for detection of *Mycoplasma* and its related strains were prepared by synthesizing 15-25 bases of specific probe with 15 bases of dT spacer at 5' end. Probes for detection of *Mycoplasma* and its related strains are not restricted to the sequences disclosed in Tables 2 and 3 and any primer
10 and probes comprising the sequences can be used in the present invention.

1. Preparation of probes for detection of *Mycoplasma* and *Ureaplasma*

15 ① Preparation of probes for genus-specific detection of *Mycoplasma* and *Ureaplasma*

For genus-specific hybridization with all *Mycoplasma* and *Ureaplasma* genus, probes of SEQ ID Nos. 7 and 8 in Table 2 were designed from conserved sequences of ITS of *Mycoplasma*. Further,
20 each Group-based conserved sequences targeted to *Mycoplasma* ITS were designed as follows. For detecting Group I (*M. arginins*, *M. arthritidis*, *M. cloacale*, *M. falconis*, *M. faucium*, *M. hominis*, *M. hyosynoviae*, *M. orale*, *M. salivarium*), probes of SEQ ID Nos. 9 and 10 were designed. For detecting Group II (*M. bovis*, *M. fermentans*, *M.*
25 *opalescens*, *M. primum*, *M. spermatophilum*, *M. synoviae*), probes of SEQ ID Nos. 11, 12, 13 and 14 were designed. For detecting Group III (*M. muris*, *M. penetrans*, *U. urealyticum*), probes of SEQ ID Nos. 15 and 16 were designed. For detecting Group IV (*M. neurolyticum*, *M. pulmonis*), probes of SEQ ID Nos. 17, 18 and 19 were designed. For
30 detecting Group V (*M. genitalium*, *M. pirum*, *M. pneumoniae*), probes of

SEQ ID Nos. 20 and 21 were designed.

② Preparation of probes for species-specific detection of *Mycoplasma* and *Ureaplasma*

For species-specific hybridization with each *Mycoplasma* and *Ureaplasma* species, 100 kind of probes of SEQ ID Nos. 28 to 127 in Table 3 were designed from species-specific sequences of ITS of *Mycoplasma* and *Ureaplasma*, which can detect 25 kind of *Mycoplasma* strains.

2. Preparation of probes for detection *Acholeplasma*

① Preparation of genus-specific probes for detection *Acholeplasma*

For genus-specific hybridization with all *Acholeplasma* genus, probes of SEQ ID No. 22 in Table 2 was designed from conserved sequences targeted to both of ITS1 and ITS2 of *Acholeplasma*. Further, each Group-based conserved sequences targeted to each *Acholeplasma* ITS1 and ITS2 were designed as follows. For Group I targeted to ITS1, probes of SEQ ID Nos. 23, 24 and 25 were designed. For Group II targeted to ITS2, probes of SEQ ID Nos. 26 and 27 were designed.

② Preparation of species-specific probes for detection *Acholeplasma*

For species-specific hybridization with each *Acholeplasma* species, probes of SEQ ID Nos. 128 to 133 in Table 3 were designed from species-specific sequences of ITS of *Acholeplasma*.

Example 3: Preparation of target DNA

1. Preparation of target DNA for detection of *Mycoplasma* and its

related strains

For preparing target DNA for detection of *Mycoplasma* and its related strains, 187~290bp size of ITS regions were selectively amplified using 5'-biotin-GTG(C/G)GG(A/C)TGGATCACCTCCT-3' (MP16SF-2) and 5'-biotin-GCATCCACCA(A/T)A(A/T)AC(C/T)CTT-3' (MP23SR-2), and 5'-biotin-AAAGTGGGCAATACCCAACGC-3' (M78) and 5'-biotin-CCACTGTGTGCCCTTTGTTTCCT-3' (R34) which were biotin-labeled respectively (Tang et al., 2000.). To prepare genomic DNAs of *Mycoplasma* and its related strains isolated in Example 1, PCR were carried out using the above primers in the following conditions: denaturation at 94°C for 3 minutes, 30 cycles of amplification at 94°C for 30 seconds, at 55°C for 2 minutes and at 72°C for 2 minutes, and final extension at 72°C for 10 minutes. After the reaction, the reaction products were analyzed by ELECTROPHORESIS on a 2% agarose gel. FIG. 3 is an electrophoresis image taken after the PCR performed using primers capable to amplify ITS target sequences of several *Mycoplasma*.

Example 4: Probe immobilization on support

Among the probes prepared in Example 2, each representative probes for *Mycoplasma*, *Acholeplasma* and *Ureaplasma* were selected. Each of the selected probes was transferred to 384-well microplate, diluted to a concentration of 50 pmole by adding spotting solution, and immobilized on a slide glass using a microarrayer (Cartesian Technologies, USA). In FIG 4, each probes for detection of *Mycoplasma* and its related strains correspond to SEQ ID Nos. 7, 28, 30, 33, 38, 41, 49, 52, 58, 61, 69, 75, 83, 85, 87, 30, 90, 92, 96, 100, 105, 110, 114, 120, 122, 22, 128, and 7 in order. Two spots of each kind of the probes were attached to the support and left in a slide box at room temperature for 24 hours or in a dry oven at 50°C for about 5 hours to be fixed to the surface of the support.

Example 5: Unimmobilized probe washing

The slide glass after the process in Example 4 was washed with a 0.2% SDS buffer solution and then distilled water at room temperature to remove unimmobilized probes. The washed slide glass was immersed in a sodium borohydride (NaBH_4) solution for 5 minutes and then washed again at 100°C . Final washing with a 0.2% SDS solution and then distilled water was followed by centrifugation to fully dry the slide glass.

Example 6: Hybridization

The biotin-labeled target products prepared in Example 3 were thermally treated to be denaturated into single strands and cooled to 4°C . A hybridization reaction solution containing $2\mu\text{l}$ of the target products was prepared. This hybridization reaction solution was portioned on the slide glass after the process in Examples 4 and 5, and the slide glass was covered with a cover slip and reacted at 25°C for 1 hours.

Example 7: Unhybridized target DNA washing

TO WASH OUT UNHYBRIDIZED TARGET DNAs, THE COVER SLIP WAS REMOVED USING A 2X SSC WASHING SOLUTION (300MM NaCl, 30MM Na-CITRATE, PH 7.0), AND THE SLIDE WAS WASHED WITH 2X SSC AND THEN 0.2X SSC, FOLLOWED BY CENTRIFUGATION TO FULLY DRY THE SLIDE GLASS.

Example 8: Staining and Result analysis

To determine hybridization of PCR products and probes, Cy5-streptavidin or Cy3-streptavidin (Amersham pharmacia biotech, USA) was diluted with 6x SSC and BSA (Bovine Serum Albumin), about $40\mu\text{l}$ of dilutes was portioned on slide glass, and the slide glass was covered

with a cover slip to block light and reacted at 50°C for about 20 minutes. After the reaction, the cover slip was removed using a 2X SSC solution, and the slide was washed with 2X SSC and then 0.2X SSC. The hybridized result was scanned using a non-confocal laser scanner (GenePix 4000A, Axon Instruments, U.S.A.) and analyzed by image analysis.

FIG. 5 shows results of image analysis of specific hybridization reaction of each probes for detecting genotypes of representative 11 kinds of *Mycoplasma* and its related strains and results of numerical analysis calculated from pixel intensity.

FIG. 5a shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 85) of *M. cloacale*. FIG. 5b shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 87) of *M. falconis*. FIG. 5c shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 90) of *M. hyosynoviae*. FIG. 5d shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 49) of *M. neurolyticum*. FIG. 5e shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 52) of *M. opalescens*. FIG. 5f shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 69) of *M. penetrans*. FIG. 5g shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 61) of *M. pirum*. FIG. 5h shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 83) of *M. salivarium*. FIG. 5i shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No.

100) of *M. spermatophilum*. FIG. 5j shows results of hybridization reaction of genus-specific probe (SEQ ID No. 7) and species-specific probe (SEQ ID No. 122) of *U. urealyticum*. FIG. 5k shows results of hybridization reaction of genus-specific probe (SEQ ID No. 22) and
5 species-specific probe (SEQ ID No. 128) of *A. laidlawii*.

Industrial Applicability

As described above, the present invention provides a rapid and accurate assay method capable of simultaneously detecting many
10 *Mycoplasma* and its related strains from a single sample using a microarray comprising novel oligonucleotides for detecting *Mycoplasma* and its related strains which are known as a source of contamination of cell lines and biological products and human pathogenic.

Also, the present invention provides an objective and credible
15 assay method capable of tracing a contamination source for preventing expansion of infective *Mycoplasma* and its related strains and controlling a contamination of *Mycoplasma* against biological products and stem cells or cord blood cells which are useful for gene therapy and cell therapy.

20 Further, the present invention provides very specific and sensitive hybridization assay for detecting *Mycoplasma* and its related strains using oligonucleotide probes designed based on sequence analysis of ITS region of many *Mycoplasma* Strains.

[References]

- 25 1. Doblhoff-Dier O, Collins CH. (2001) Biosafety: future priorities for research in health care. J Biotechnol. 85 : 227-39.
2. Jung H, Wang SY, Yang IW, Hsueh DW, Yang WJ, Wang TH, Wang HS. (2003) Detection and treatment of *Mycoplasma* contamination in
30 cultured cells. Chang Gung Med J. 26 : 250-8.
3. Wisher M. (2002) Biosafety and product release testing issues

- relevant to replication-competent oncolytic viruses, Review. Cancer Gene Ther. 9 : 1056-61.
4. Hopert A, Uphoff CC, Wirth H, Hauser H, Drexler HG. (1993) Specificity and sensitivity of polymerase chain reaction in comparison with other methods for the detection of *Mycoplasma* contamination in cell lines. J Immunol Methods. 164 : 91-100.
 5. Kong F, James G, Gordon S, Zelyski A, Gilbert GL. (2001) Species-Specific PCR for Identification of Common Contaminant Mollicutes in Cell Culture. Appl Environ Microbiol. 67 : 3195-200.
 - 10 6. Dorigo-zetsma J.W., Zaat SAJ, Wertheim-van D, Spanjaard PME, Rijntjes J, Waveren V, Jensen JS, Angulo AF, Dankert J. (1997) Comparison of PCR, culture, and serological tests for diagnosis of *Mycoplasma pneumoniae* respiratory tract infection in children. J Clin Microbiol. 37 : 14-7.
 - 15 7. Jensen JS, Borre MB, Dohn B. (2003) Detection of *Mycoplasma genitalium* by PCR Amplification of the 16S rRNA Gene. J Clin Microbiol. 41 : 261-266.
 8. Uphoff CC, Drexler HG. (2002) Comparative PCR analysis for detection of *Mycoplasma* infections in continuous cell lines. In Vitro Cell Dev Anim. 38 : 79-85.
 - 20 9. Gohlman HW, Weiner J 3rd, Schon A, Herrmann R. (2000) Identification of a small RNA within the pdh gene cluster of *Mycoplasma pneumoniae* and *Mycoplasma genitalium*. J Bacteriol. 182 : 3281-4.
 - 25 10. Tang J, Hu M, Lee S, Roblin R. (2000) A polymerase chain reaction based method for detection of *Mycoplasma/Acholeplasma* contaminants in cell culture. J Microbiol Methods. 39 : 121-6.

What is claimed is:

1. An ITS (Internal transcribed spacer) target DNA for genotyping *Mycoplasma* strains, comprising any one sequence selected from SEQ ID Nos. 1 to 6.
5
2. An oligonucleotide for genus-specific genotyping of *Mycoplasma* and *Ureaplasma* strains, comprising any one sequence selected from SEQ ID Nos. 7 to 21 or its complementary sequence.
10
3. An oligonucleotide for species-specific genotyping of *Mycoplasma* and *Ureaplasma* strains, comprising any one sequence selected from SEQ ID Nos. 28 to 127 or its complementary sequence.
- 15 4. An oligonucleotide for genus-specific genotyping of *Acholeplasma* strains, comprising any one sequence selected from SEQ ID Nos. 22 to 27 or its complementary sequence.
- 20 5. An oligonucleotide for species-specific genotyping of *Acholeplasma* strains, comprising any one sequence selected from SEQ ID Nos. 128 to 133 or its complementary sequence.
- 25 6. A microarray comprising more than one oligonucleotide selected from genus-specific and species-specific oligonucleotides for genotyping *Mycoplasma*, *Acholeplasma* and *Ureaplasma* strains according to any one from claims 2 to 5 as probes attached on a support.
- 30 7. The microarray according to claim 6, wherein the probes are any one selected from a group consisting of DNA, RNA, PNA, LNA and HNA.
8. The microarray according to claim 6, wherein the support is any

one selected from a group consisting of slide glass, plastic, membrane, semiconductive chip, silicon and gel.

9. A method for detecting *Mycoplasma* strains, comprising the
5 following steps:
- a) extracting nucleic acids from a sample;
 - b) amplifying target DNA among the extracted nucleic acids;
 - c) hybridizing the amplified target DNA with probes of the
microarray according to claim 6; and
 - 10 d) detecting signals generated from the hybridization reaction.

10. A kit for diagnosing *Mycoplasma* infection, comprising more than
one oligonucleotide selected from genus-specific and species-specific
oligonucleotides for genotyping *Acholeplasma*, *Mycoplasma* and
15 *Ureaplasma* strains according to any one from claims 2 to 5.

Abstract of the Invention

The present invention relates to a method for detecting *Mycoplasma* and its related strains which are source of contamination of cell lines and biological products and human pathogenic. More particularly, the present invention relates to genus-specific and species-specific oligonucleotides for genotyping of *Mycoplasma*, *Acholeplasm* and *Ureaplasma* strains, microarray comprising the oligonucleotides, and method for detection of species using the microarray.

As described above, the present invention provides a rapid and accurate assay method capable of simultaneously detecting many *Mycoplasma* and its related strains from a single sample using a microarray comprising novel oligonucleotides for detecting *Mycoplasma* and its related strains which are known as a source of contamination of cell lines and biological products and human pathogenic. Further, the present invention provides an objective and credible assay method capable of tracing a contamination source for preventing expansion of infective *Mycoplasma* and its related strains and controlling a contamination of *Mycoplasma* against biological products and stem cells or cord blood cells which are useful for gene therapy and cell therapy.

20

FIG. 1a

<i>M. bovis</i>	-----ATA-----	TGTTCTTTGAAAACTGAATAGTAAATATTTTT	142
<i>M. primetus</i>	-----TT-----	TGTTCTTTGAAAACTGAATAGTAAATATTTTT	181
<i>M. farinifera</i>	-----ATT-----	TGTTCTTTGAAAACTGAATAGTAAA---TTTTT	177
<i>M. upatseuensis</i>	-----TT-----	TGTTCTTTGAAAACTGAATAGTAAA---TTTTA	159
<i>M. apicostrophilus</i>	-----AT-----	TGTTCTTTGAAAACTGAATAGTAAA---TTTTT	196
<i>M. aynoviae</i>	-----AAI-----	TGTTCTTTGAAAACTGAATAGTAAA---TAA	100
<i>M. neurolyticus</i>	TAATAAATGTTTT-----AAT	TGTTCTTTGAAAACTGAATAGTAAA---TA--T	176
<i>M. pulchra</i>	---AACAAATA-----	TGTTCTTTGAAAACTGAATAGTAAA---TAAAT	159
<i>M. hyorhiniae</i>	-----ATA-----	TGTTCTTTGAAAACTGAATAGTAAA---TAA	112
<i>M. wiluifidus</i>	-----TT-----	TGTTCTTTGAAAACTGAATA-----T	115
<i>M. faucae</i>	-----TT-----	TGTTCTTTGAAAACTGAATA-----T	123
<i>M. orele</i>	-----II-----	TGTTCTTTGAAAACTGAATA-----I	108
<i>M. hyoaynoviae</i>	-----TT-----	TGTTCTTTGAAAACTGAAT-----T	119
<i>M. salivarius</i>	-----TT-----	TGTTCTTTGAAAACTGAAT-----T	115
<i>M. falconia</i>	-----TT-----	TGTTCTTTGAAAACTGAATA-----T	92
<i>M. fusaria</i>	-----	TGTTCTTTGAAAACTGAATA-----T	87
<i>M. arginini</i>	-----TT-----	TGTTCTTTGAAAACTGAATA-----T	93
<i>M. oloaeae</i>	-----II-----	TGTTCTTTGAAAACTGAATA-----I	96
<i>M. granitarius</i>	CCAGTTCTGAAAG--AATGTTTTTGAA	TGTTCTTTGAAAACTGAATA-----T	160
<i>M. pneumoniae</i>	CCAGTTCTGAAAG--AACATTTCCGC	TGTTCTTTGAAAACTGAATA-----T	190
<i>M. pirus</i>	TAAATTTTAAABTAGTAGAGATGG	TGTTCTTTGAAAACTGAATA-----T	213
<i>M. muris</i>	TT-----	CTTTGAAAACTGAATATTTGTA-----	106
<i>M. panatrans</i>	TT-----	CTTTGAAAACTGAATATTTATA-----	184
<i>U. urealyticus</i>	IIAAIIIAIATG--GATGATGGA	CTTTGAAAACTGAATATTTATA-----	199

**** * * * * *

FIG. 1b

<i>M. arthritidis</i>	AAA—CATCGTATCCAGTTTGGAGAACTAAACTTCCTCTCTTTGTTCTTTGAAAC	108
<i>M. taucium</i>	AAA—CATCGTATCCAGTTTGGAGAACTAAACTTCCTCTCTTTGTTCTTTGAAAC	88
<i>M. falconis</i>	TAA—TATTGTATCCAGTTTGGAGAACTA—ACTCTCTTTT—GTTCTTTGAAAC	85
<i>M. hominis</i>	AAAAAATTATTGTATCCAGTTTGGAGAACTTA—TCTCTCG—GTTCTTTGAAAC	80
<i>M. arginini</i>	AAA—TATTGTATCCAGTTTGGAGAACTA—TCTCTCAATTT—GTTCTTTGAAAC	86
<i>M. cloacale</i>	GAATTAATTGTATCCAGTTTGGAGAACTT—TCTCTCAATTTGTTCTTTGAAAC	89
<i>M. hyosynoviae</i>	CA—TATCGTATCCAGTTTGGAGAACTTAT—TCTCTCTTTT—GTTCTTTGAAAC	113
<i>M. orale</i>	CAA—TATCGTATCCAGTTTGGAGAACTAT—CTCTCAATTT—ATTCTTTGAAAC	102
* * ***** *		
<i>M. arthritidis</i>	—TTAAAAAATTAATATTTCAAA—GTTTAGATCAAOCATAGAAATACAA	173
<i>M. taucium</i>	—TTAAAAAATTAATATTTCAAA—GTTTAGATCAAOCATAGAAATACAA	153
<i>M. falconis</i>	ATTA—TTAATTAATATTTCAAA—GTTTAGATCAAOCATAGAAATACAA	150
<i>M. hominis</i>	—TA—TTAATTAATATTTCAAA—GTTTAGATCAAOC—ATAGAAATATTT	141
<i>M. arginini</i>	ATTAAATTTATTAATATTTCAAA—GTTTAGATCAAOCATAGAAATATAT	153
<i>M. cloacale</i>	—TCAATAAATTAATATTTCAAA—GTTTAGATCAAOCATAGAAATATTC	154
<i>M. hyosynoviae</i>	A—TTATCAAATTAATATTTCAAA—GTTTAGATCAAOCATAGAAATATTC	178
<i>M. orale</i>	—TTAAAAATTAATATTTCAAAA—ATTTAGATCAAOCATAGAAATATTC	166
***** * ***** * *****		

2/10

FIG. 1c

<i>M. bovis</i>	TTATTAAATAGGTCAAAAGCTA	ATATCTAGTTTTGAAGAGCA	TTCTCTCAT	144
<i>M. primatum</i>	TT--TAATAGGTGAGGAGCTT	ATATCTAGTTTTGAAGAGCA	TTCTCTCTT	148
<i>M. fermentans</i>	TTTTTATGAGGTCTAAAGCTT	ATATCTAGTTTTGAAGAGCA	ATATTTTCTCTCAT	148
<i>M. opalescens</i>	T-----ATGATCTACAAAGCT	ATATCTAGTTTTGAAGAGCA	TTCTCTCTT	129
<i>M. spermatophilum</i>	TT--TTATGAGGTCTAAAGCTT	ATATCTAGTTTTGAAGAGCA	TCTCTCTAAT	166
<i>M. synoviae</i>	GCTTTTTTGGCTTGGCTAT	GTATTTAGTTTTGAAGAGCA	CTCTCTTAAA	141

* * * *** ***** ** ***

<i>M. bovis</i>	ATGTTCTTTGAAAAGTGAATAGTAAATATTTTTC	GATATTTACAGGACATGAAA	201
<i>M. primatum</i>	-TGTTCCTTTGAAAAGTGAATAGTAAATATTTTTC	GATATTTACAGGACATGCACTC	207
<i>M. fermentans</i>	TGTTCCTTTGAAAAGTGAATAGTAA--TTT	GATATTTACAGGACATGAAA	200
<i>M. opalescens</i>	-TGTTCCTTTGAAAAGTGAATAGTAA--TTT	GAATATTTACAGGACATGATA	182
<i>M. spermatophilum</i>	-TGTTCCTTTGAAAAGTGAATAGTAA--TTT	TTTATTTACAGGACATGATA	219
<i>M. synoviae</i>	TGTTCCTTTGAAAAGTGAATAGTAA--TAA	GATATTTACAGGACATGAAAAT	193

***** ***** * ***** ** * * * *

<i>M. bovis</i>	ATCAA--TTAA	GGTGAATTGTTTGGAT	CATCGAGT	AAGTCATATTTA	250
<i>M. primatum</i>	CCATCAAAAATTTAA	GGGGAATTGTTTGGAT	CATCGAGA	AAATCATATTAA	261
<i>M. fermentans</i>	--TTAA--TTAA	GGTGAATTGTTTGGAT	TCATCGAGA	AAATCATATTAA	250
<i>M. opalescens</i>	ATTAAATGATTTA	GGTGAATTGTTTGGAT	CATCGAGATA	AAAAACAATCATAA	238
<i>M. spermatophilum</i>	TAATTGAA--TTAA	GGTGAATTGTTTGGAT	TCATCGAGA	TAGTCATTTTAA	270
<i>M. synoviae</i>	ATAAATTAATTAA	GGTGAATTGTTTGGAT	ACCGAGTT	TAAATTAT-TGAA	243

* * * * * ***** *** * **** * ** *

<i>M. bovis</i>	TATGATTCATGGAATGCTT	AAATACACATCTAA	ACTAACCAATAGGA	304
<i>M. primatum</i>	TATGATTCATGGAATGCTT	AAATACACATCTTAA	ACTAA--ACAATAGGA	319
<i>M. fermentans</i>	TATGATTCATGGAATGCTT	AAATACACATCATAACA	AACATAACCAATAGGA	306
<i>M. opalescens</i>	TTTTTGATTCATGGAATGCTT	AAATACACATCATAATGT	AACCAATACCAATAGGA	296
<i>M. spermatophilum</i>	AATGATTCATGGAATGCTT	AAATACACATCAAAACAA	ACAATCTATACCAATAGGA	330
<i>M. synoviae</i>	AATATTTATTAATGCTT	GAATACA--TCATAAC	AATATAACCAATAGGA	295

* *** *** ***** ***** ** ** *****

FIG. 1d

```

M. muris      CCTCCTTTCTATCGGAGTACA TTTAGATTATTACACCATATTAGAATATTTTAAATATT  50
M. penetrans CCTCCTTTCTATCGGAGTACA TTAAGCTAAGTAACAAATATTAG-----ATATATT  52
U. urealyticum CCTCCTTTCT-TGGAGTAAAT TTTTAAT--TTACGTACTAATAAG-----TGTACATTTT  53
*****

```

```

M. muris      TGTGTACTTT-TTATAGAAAACCCCCACATCAATAAAGCTAA-----ATAAAAAATTATT  115
M. penetrans TGTGTACTTTATTAATAAAAAATCCCTAAAGTGAATTTATCTCATGTTATATAAGAGTAAGT  112
U. urealyticum TATTAAAAATCCATATGAATATAAGCCACTTTTTTAAAAATTTT-----TAAAAATTCATAT  109
* * * * *

```

```

M. muris      TTTCGG-CGGATTCTATTTAGTTTTCAGAGATA-TTTCTCTCATGATAGTT-----  165
M. penetrans TCTAGG-CGGATTCTATTTAGTTTTCAGAGAT-TTTCTCTCAAGATAGTT-----  162
U. urealyticum ---GG-CGGATTCTATTTAGTTTTCAGAGTTTATTCTCTCCCATATAATTAAATTTATT  165
** *****

```

3/10

FIG. 1e

M. pulmonis CTACGGAGTACAAAACCATTTTTTTAGAAATGGCATTYY-----TCTATCAATAGTTAT- 54
M. neurolyticum CTACGGAGTACACATACATCTTATTAATTTGGTTATTTAAAAATCCTTTTATAATAAAT 60
 ***** * *** ** *** + **** ** ** * *** *

M. pulmonis -AGAAAGTCCTTATGTGTATCTTCCCAATTAGATATCTAGTATTCACTTTTGAAAGTCT 113
M. neurolyticum AAAAAGGTTATTATGGG-CTTCCCAATAG-TTTTGTATCTAGTTTGAAAGATTT 114
 * * * * * ***** * ***** * * * * * ***** * *

M. pulmonis A-----TCTTTCAAA-ACAAATA-----GTCTTTAAAAACTGAATAGCATAT 155
M. neurolyticum AATTTTTTCTTTCTAATTAATAAATGTTTTAATATATTCTTTGAAAACTGAATAGCAAT 174
 * ***** ** * **** ***** ***** *

M. pulmonis AAATTAATATGATAACGTCATCAAAATGTAAATTTTTTGATCCGAGTCATTTTTTAACAA 215
M. neurolyticum ---ATTGAAATTTTAATTTTCATAATATTTCAACAACGACATTACAACAACGAGTCTAACTG 232
 *** * ** *** ***** * * * * * * * * * * * *

M. pulmonis TTT-GTTAAAAAAT-AAAATAGTACCTTAGG-ATAACATCAAAAA-ATAAAT 265
M. neurolyticum TTTTATTGAACAGTTAGCTTAAATAGTACCTTAGATATATAAATCTAAAACAATAGGC 292
 *** ** * * * ***** * * * * * * * * * * *

FIG. 1f

M. pneumoniae AACATTTCCGC-----TTCTTTCAAAACTGAAACGACAA-TCTTTCTAGTTCA----- 205
M. genitalium AATGTTTTTGAACAGTTCTTTCAAAACTGAAACGACAA-TCTTTCTAGTTCA----- 175
M. pirum AGTAGAGATGG-----TTCTTTGAAAACTGGAATACACAAATCTTTCTAGTTCTTTGTGTG 235
 * * ***** * * * * * *****

M. pneumoniae ---AA-TAAATACCAAAGG-ATCAATAC-AATAAGTTACTAAGGGCTTATGCT 252
M. genitalium ---AAAATAAATACCAAAGG-ATCAATAC-AATAAGTTACTAAGGGCTTATGCT 224
M. pirum AATAACACAAATATCAAATGCTAATGGATATCAAAATAAGTTACTAAGGGCTTATGCT 295
 ** ***** ***** * ** *** *****

[illegible]

<i>A. tridivarii</i>	1	CAAGTAACCAATATTATAATAAGT	GGGGCTGTAGCTCAGTTGTTAGAGGCTGTGCT	168
<i>A. denii</i>	1	CAA-----AAGT	GGGGCTGTAGCTCAGTTGTTAGAGGACAGCT	155
<i>A. exanthus</i>	1	TAG-----TAA	GGGGCTGTAGCTCAGTTGTTAGAGGACAGCT	133
<i>A. rodjensis</i>	1	TTA-----	TGGGCTATAGCTCAGCTGTTAGAGGACAGCC	122

```

A. tridivii 1      TGATAAAGGTGAGGTGATGGTTCAGTCCGTTGAGGCCACCATTAATAAATATCAATA 227
A. oculi 1         TGATAAAGGTGAGGTGATGGTTCAGTCCATTCAGGCCACCAT----- 201
A. tridivii 1     TGATAAAGGTGAGGTGATGGTTCAGTCCATTCAGGCCACCATTTATAT----- 184
A. modicum 1      TGATAAAGGTGAGGTGATGGTTCAGTCCATTCAGGCCACCATTATAG----- 172

```

A. f. id fawii GTAA—TATTCTCTAAATTGTCATCATATTCACTTTTGAAAGACTTAA—AGTAATF— 104
A. oculi GCAAGCAATTCTCTAA—TTTGTTCATCATATTCACTTTTGAAAGACTTAA—TCCAAGTG 115
A. exanthema —AAAGCAATTTCTTCA—TTTGTTCATCATATTCACTTTTGAAAGACTTTG—ACTTGTFF 98
A. modicus —————CATTTTCATCATATTCACTTTTGAAAGATTTTCTTCTAATAT 84
 * ***** ** **

[illegible]

FIG. 3a

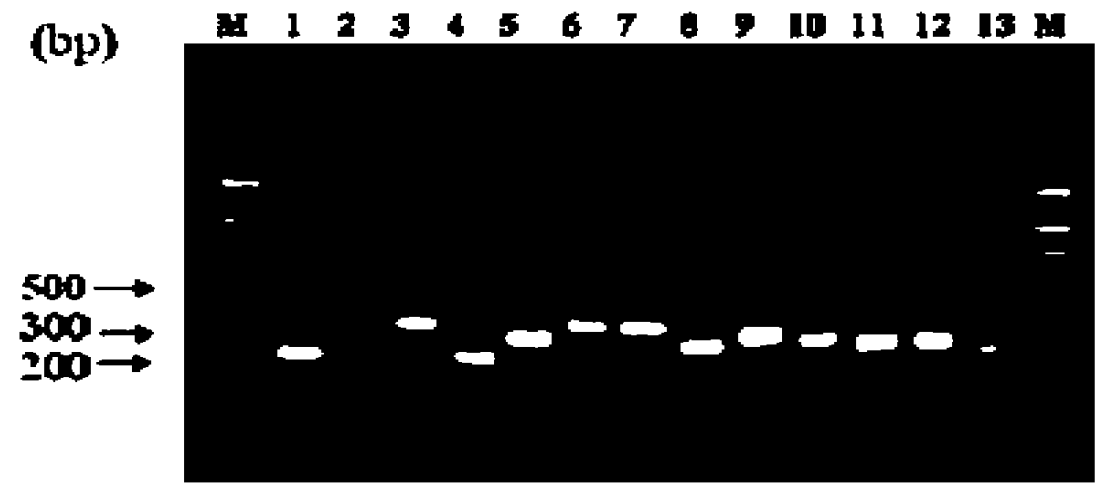


FIG. 3b




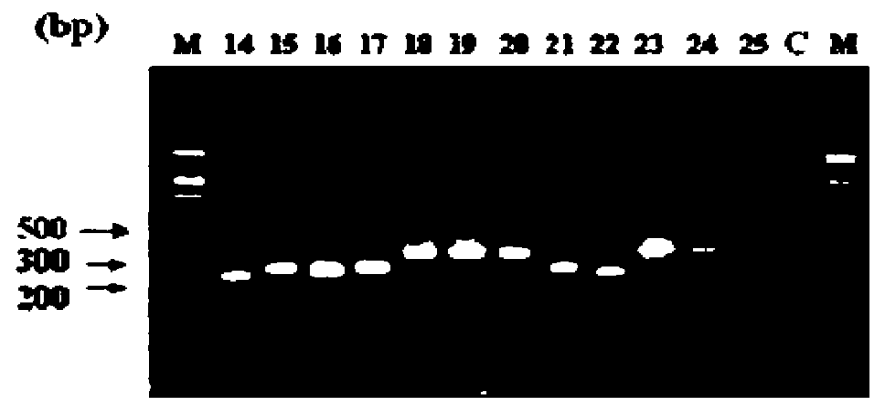



FIG. 4


MP-C [7]	<i>M. arginini</i> [28]	<i>M. arthrodia</i> [30]	<i>M. fermentans</i> [33]	<i>M. hominis</i> [38]	<i>M. hyorhinis</i> [41]
<i>M. neurolyticum</i> [49]	<i>M. opalescens</i> [52]	<i>M. orale</i> [58]	<i>M. pirum</i> [61]	<i>M. penetrans</i> [69]	<i>M. pulmonis</i> [75]
<i>M. salivarium</i> [83]	<i>M. cibacale</i> [85]	<i>M. faconis</i> [87]	<i>M. faucium</i> [30]	<i>M. hyosynoviae</i> [90]	<i>M. muris</i> [92]
<i>M. primatum</i> [96]	<i>M. spermatophilum</i> [100]	<i>M. synoviae</i> [105]	<i>M. pneumoniae</i> [110]	<i>M. genitalium</i> [114]	<i>M. bovis</i> [120]
<i>U. urealyticum</i> [122]			AP-C [22]	<i>A. laidlawii</i> [128]	MP-C [7]

*[] corresponds to SEQ ID No's of Tables 2 and 3.

[illegible]

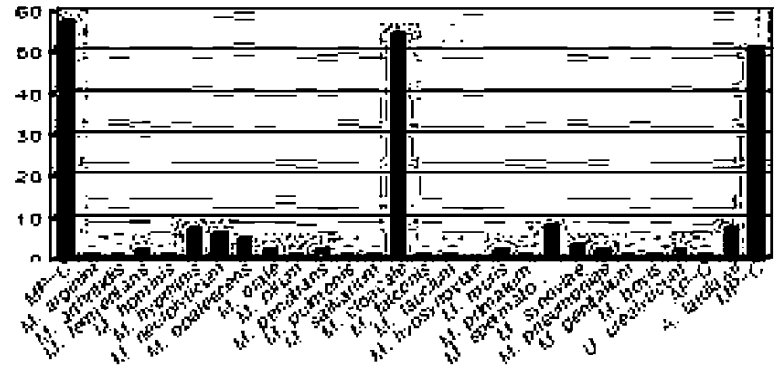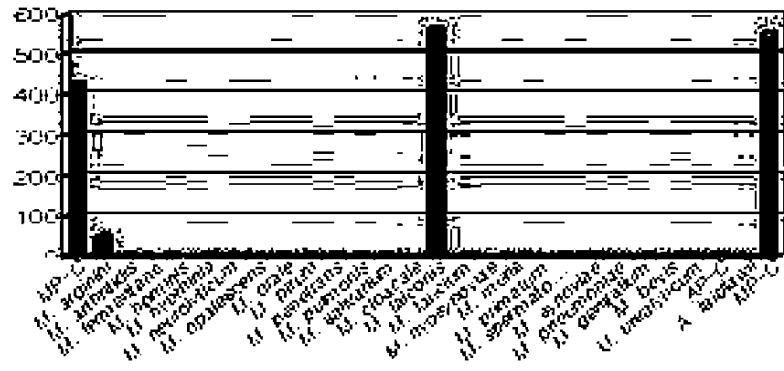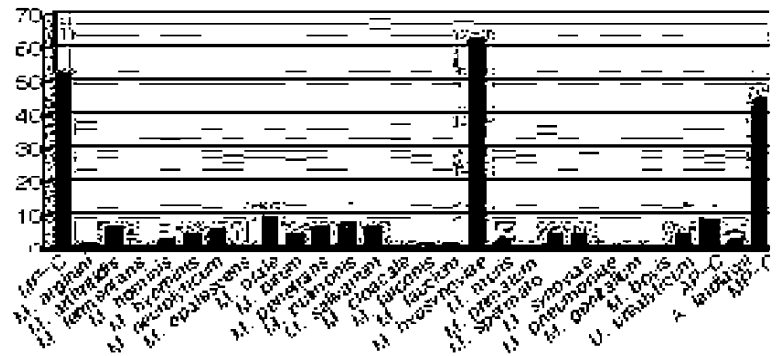

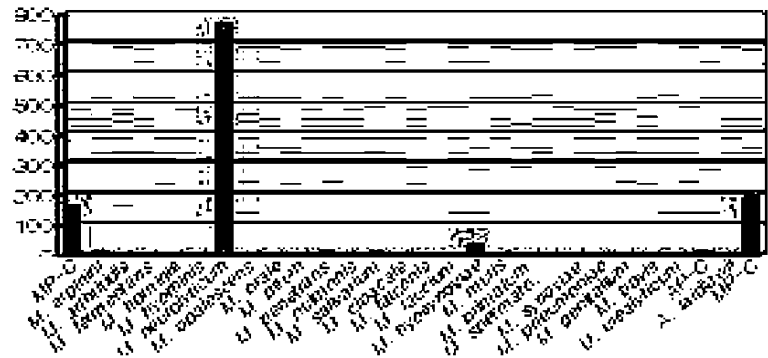

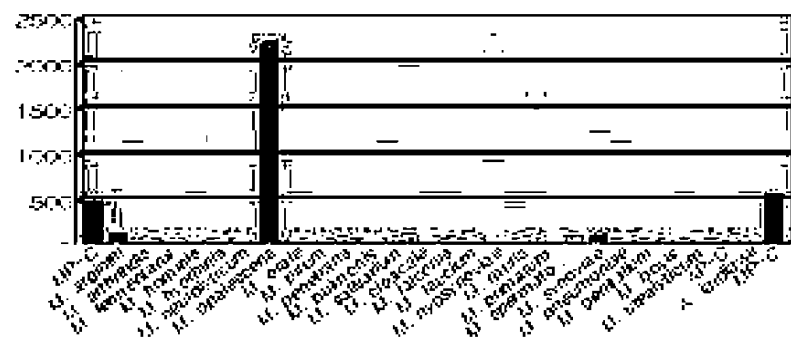

the 1990s, the number of people in the United States who are aged 65 and older has increased by 25% (U.S. Census Bureau, 1997). The number of people aged 65 and older is projected to increase by 50% by the year 2020 (U.S. Census Bureau, 1997). The number of people aged 65 and older is projected to increase by 50% by the year 2020 (U.S. Census Bureau, 1997).

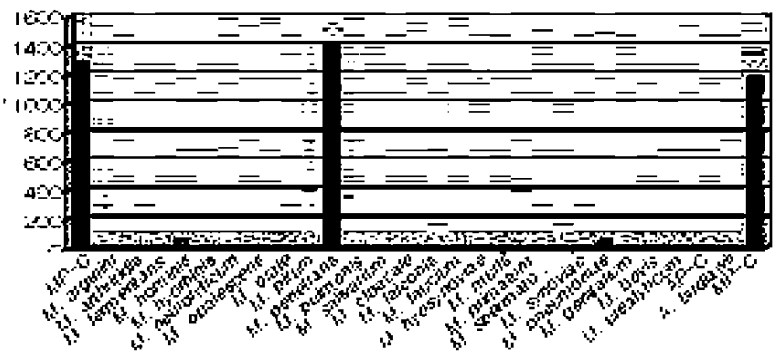


FIG. 5g

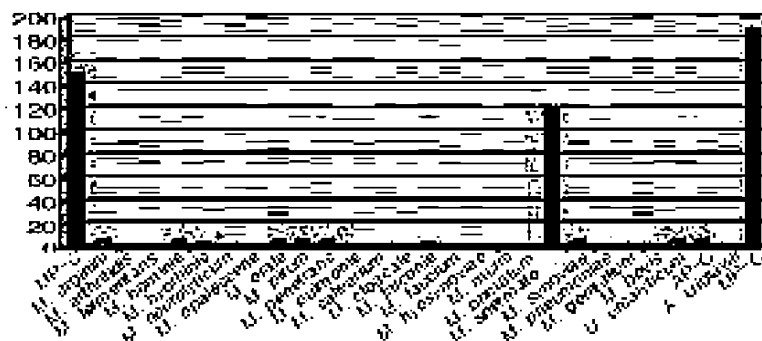
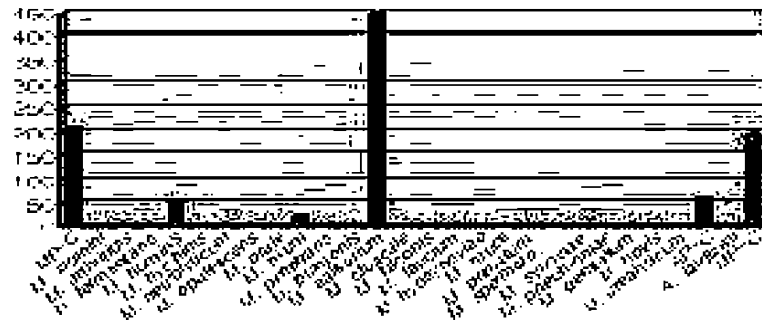

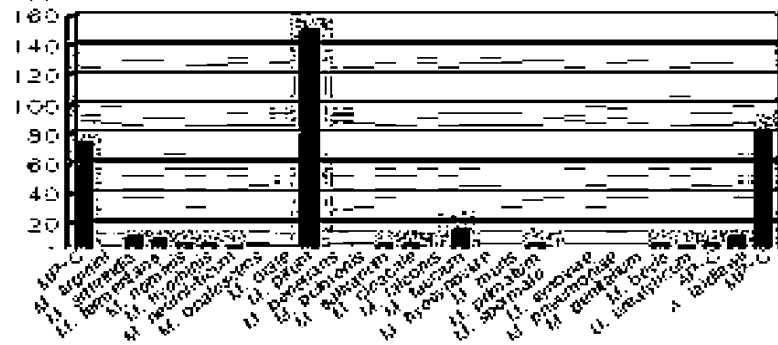


FIG. 5j

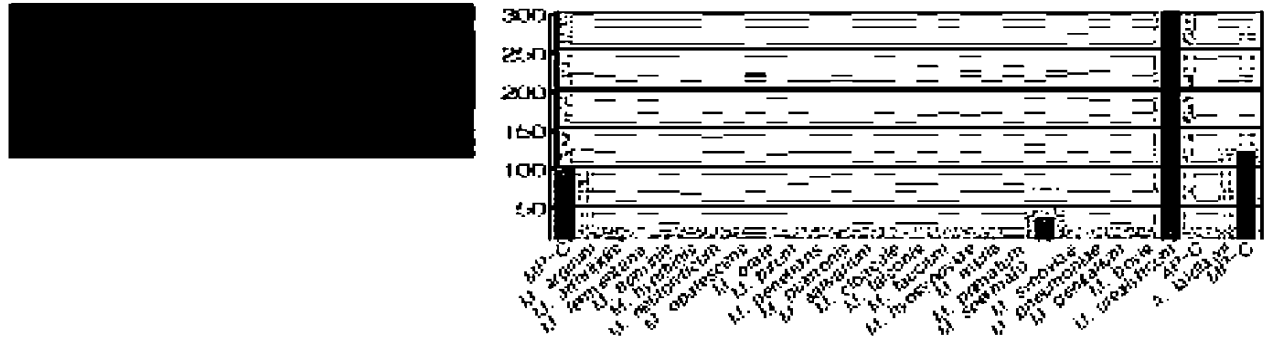


FIG. 5k

